

ToeUP™ - a lightweight ankle foot orthosis for Foot Drop

A dynamic dorsiflexion assist orthosis from Allard

ToeUP™ is an ankle foot orthosis designed for individuals with mild foot drop. The soft textile offers comfort around the ankle and combined with the integrated elastic cords provides stability and dynamic support to help lift the forefoot during walking.

- The support is adjusted by the placement of hook and elastic cord
- Soft and flexible
- For outdoor and indoor use
- The D-ring can hold the elastic bands together for a more discrete profile



Dynamic dorsiflexion assist, support, and comfort for mild foot drop

ToeUP™ features a low-profile component that wraps around the ankle just above the malleoli. Elastic cords with hooks connect the ankle section to the shoe, secured either beneath the laces or through the lacing holes. This system helps lift the forefoot during ambulation. The amount of lift can be easily adjusted by changing the position of the hooks or the strap on the ankle section.

The supplementary D-ring can be used to hold the elastic cord together in desired location.

ToeUP™ includes a forefoot wrap that may be used to attach the elastic cord when not wearing shoes or unlaced shoes, such as sandals.



The D-ring holds the straps together.




Forefoot wrap for indoor use or with unlaced shoes.



Item no.	Model	Size	L/R	Circumference of the ankle (1)	Circumference around foot (2)
285000010	ToeUP	X-Small	Bilat.	13-17 cm	-
285000011	ToeUP	Small	Bilat.	17-21 cm	-
285000012	ToeUP	Medium	Bilat.	21-25 cm	-
285000013	ToeUP	Large	Bilat.	25-29 cm	-
285011011	ToeUP Foot Part	Small	Left	-	17-21 cm
285011012	ToeUP Foot Part	Medium	Left	-	21-25 cm
285011013	ToeUP Foot Part	Large	Left	-	25-29 cm
285012011	ToeUP Foot Part	Small	Right	-	17-21 cm
285012012	ToeUP Foot Part	Medium	Right	-	21-25 cm
285012013	ToeUP Foot Part	Large	Right	-	25-29 cm
285020011	ToeUP Elastic Cord	X-Small - Small	-	-	-
285020013	ToeUP Elastic Cord	Medium - Large	-	-	-
285030000	ToeUP Hooks	-	-	-	-
285040000	ToeUP D-ring	-	-	-	-
285050000	Under Lace Hook	-	-	-	-
285060010	ToeUP Ankle Part	X-Small	Bilat.	13-17 cm	-
285060011	ToeUP Ankle Part	Small	Bilat.	17-21 cm	-
285060012	ToeUP Ankle Part	Medium	Bilat.	21-25 cm	-
285060013	ToeUP Ankle Part	Large	Bilat.	25-29 cm	-

Please note that item 28500 includes the following components: 28502, 28503, 28504, 28505, and 28506. These parts are included when ordering the main product (28500). Items 28501–28506 can also be ordered separately as spare parts if needed.

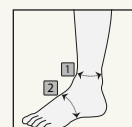



Support for Better Life

Everyone should be able to live their life to the fullest, regardless of their mobility challenges. With innovative solutions developed in close collaboration with healthcare professionals and patients, we strive to provide Support for Better Life.

allard

Allard USA

300 Forge Way, Suite 3, Rockaway, NJ 07866-2056, USA | +1 888 678 6548 | info@allardusa.com | allardusa.com