

## S.O.T THUMB ORTHOSIS

Positions the thumb in an optimal resting position.



Soft yet compressing textile cover with no hard seams or edging.



Molded and shaped palmar section with aluminum reinforcement to support the thumb in an optimal resting position.

- Designed to restrict movement in the CMC and MCP joints
- Support the hand's arch
- Position the thumb's joints and muscles in a comfortable/desirable resting position
- Avoid overstretched ligaments that can lead to joint subluxation



# S.O.T THUMB ORTHOSIS


is a completely new anatomically shaped orthosis for thumb and palm.

## INDICATIONS

Osteoarthritis, rheumatoid arthritis, adduction contracture of the thumb, overload injury of the CMC joint.




## LIMIT MOVEMENT

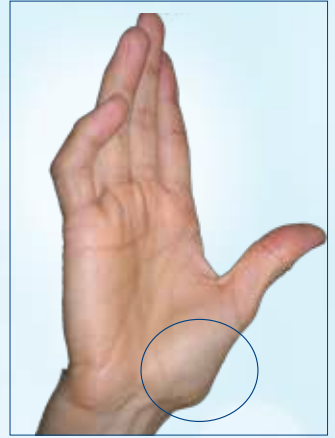

- Adjustable aluminum core can easily be customized/formed to fit individual anatomy
- Can be adjusted as the contracture changes
- Place in laundry bag, with cover attached, and machine wash at 104°F
- Excellent alternative to custom made thumb orthosis

Item No. Left	Item No. Right	Description	Size	MCP Width**	Length * * (from wrist to fingertip)
28720 1011	28720 2011	S.O.T Thumb Orthosis*	Small	≤3"	≤ 7.25"
28720 1012	28720 2012	S.O.T Thumb Orthosis*	Medium	≤ 3.25"	≤ 7.75"
28720 1013	28720 2013	S.O.T Thumb Orthosis*	Large	≤ 3.5"	≤ 8.25"
28721 1011	28721 2011	Extra cover	Small		
28721 1012	28721 2012	Extra cover	Medium		
28721 1013	28721 2013	Extra cover	Large		

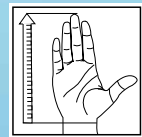

\* complete orthosis with one cover  
 \*\* Relation between MCP width and wrist to fingertip length can vary. Regard the size interval only as a guide.